




Save the Gonads shields are specifically designed to protect the reproductive organs of babies receiving X-rays.

Blocks harmful radiation

- Shields are composed of 1 mm thick lead to protect babies' gonads from the harmful radiation associated with X-rays
- Helps clinicians comply with JCAHO and FDA patient radiation safety standards¹


Designed for comfort and flexibility

- Small, lightweight design and soft satin cover ensure comfort for baby's fragile skin
- Unique heart shape allows placement on both males and females
 - Oriented upward for girls and downward for boys

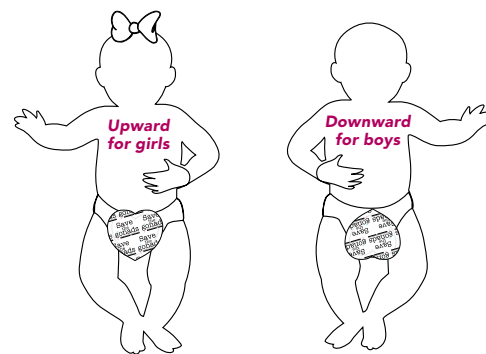


Meets infection control standards

- Designed for single patient use
- Shields are fully encased in satin to facilitate cleaning with alcohol between X-rays
 - Can be used repeatedly on the same baby

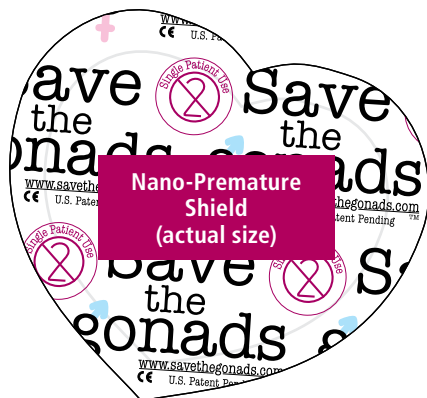

JCAHO and Federal regulations (US FDA Code of Federal Regulations Title 21, Chapter 1, Sec. 1000.50) recommend gonadal shielding for all patients – and many states require it.

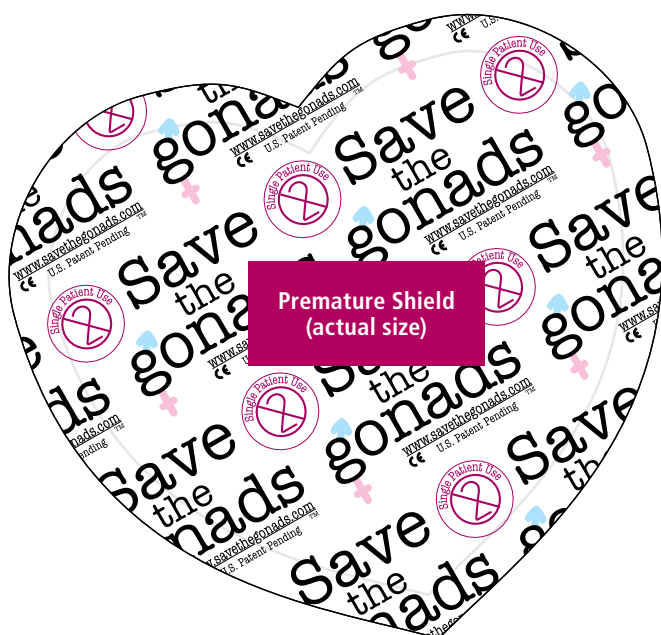

Premature babies receive an average of 30 X-rays during their NICU stay.¹ Save the Gonads shields provide added safety and protection needed for a healthy development.

Economically priced and packaged in 4 sizes

| Size | Part number | Packaging |
|-----------------|-------------|-----------------------|
| Nano-Premature | 900835 | 25 units/pack |
| Micro-Premature | 900772 | 25 units/pack |
| Premature | 900773 | 25 units/pack |
| Newborn | 900774 | 10 units/pack |
| Sample Pack* | 900775 | 10 units/variety pack |

*Sample pack includes 3 Nano-Premature, 3 Micro-Premature, 3 Premature, and 1 Newborn


Technical specifications

| | | Nano-Premature | Micro-Premature | Premature | Newborn |
|--------|-----|----------------|-----------------|-----------|---------|
| A | mm | 45.7 | 55.9 | 76.2 | 94.0 |
| | in. | 1.8 | 2.2 | 3.0 | 3.7 |
| B | mm | 55.9 | 66.0 | 86.4 | 104.1 |
| | in. | 2.2 | 2.6 | 3.4 | 4.1 |
| C | mm | 39.4 | 50.8 | 68.6 | 83.8 |
| | in. | 1.55 | 2.0 | 2.7 | 3.3 |
| D | mm | 49.5 | 61.0 | 86.4 | 94.0 |
| | in. | 1.95 | 2.4 | 3.4 | 3.7 |
| Lead | mm | 1.07 | 1.07 | 1.07 | 1.07 |
| | in. | .042 | .042 | .042 | .042 |
| Total | mm | 1.35 | 1.35 | 1.35 | 1.35 |
| | in. | .053 | .053 | .053 | .053 |
| Weight | gm | 16 | 24.6 | 47.3 | 71.7 |
| | oz. | 0.60 | 0.87 | 1.67 | 2.53 |



"Using STGs is an effective way to demonstrate your commitment to exceptional care and the attention to detail that you provide for your patients."

Pam Krieg, RN
Creator, Save the Gonads™

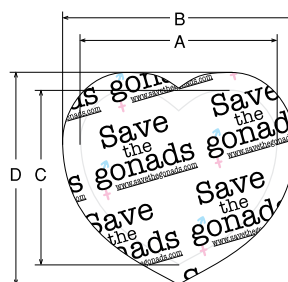

Natus Medical Incorporated
1501 Industrial Road
San Carlos, CA 94070 USA
1-800-303-0306
+1-650-802-0400
www.natus.com

1 Reference on file at Natus
Save the Gonads is a trademark of Save the Gonads Ltd. U.S. Patent Pending. www.savethegonads.com

© 2013 Natus Medical Incorporated. All Rights Reserved. All product names appearing on this document are trademarks or registered trademarks owned, licensed to, promoted or distributed by Natus Medical Incorporated, its subsidiaries or affiliates.