

SpineJet XL Fusion Handpiece Specifications

Description	Catalog Number	Procedure(s)	Working Shaft Length	Cutting Edge Width	Cutting Edge Height	Shaft Type	Тір Туре
20° XLS	57662	Open TLIF/PLIF/ALIF	20 cm	7 mm	4 mm	Straight	Rectangular
75° XLS	57667	Open TLIF	18.4 cm	7 mm	4 mm	Straight	Rectangular
20° XLS	57162	Open TLIF/PLIF/ALIF	20 cm	7 mm	4 mm	Straight	Oval
75° XLS	57167	Open TLIF	18.4 cm	7 mm	4 mm	Straight	Oval
20° XL MIS	58662	MISTLIF/PLIF	19.5 cm	7 mm	4 mm	Bayonet	Rectangular
30° XL MIS	58663	MISTLIF	19.4 cm	7 mm	4 mm	Bayonet	Rectangular
20° XL MIS	58162	MISTLIF/PLIF	19.5 cm	7 mm	4 mm	Bayonet	Oval
30° XL MIS	58163	MISTLIF	19.4 cm	7 mm	4 mm	Bayonet	Oval
20° XLE	59662	DLIF	22.5 cm	7 mm	4 mm	Straight	Rectangular
10° DL (pull)	59961	DLIF/ALIF	22.5 cm	9 mm	5 mm	Straight	Oval
10° DL (push)	59971	DLIF	22.5 cm	9 mm	5 mm	Straight	Oval

All products are disposable, single use, supplied sterile.

SpineJet XL Fusion Accessories

Description	Catalog Number		
SpineJet XL Quick Connect* (Single use, supplied sterile, for use with all Blue handle handpieces)	55500		
SpineJet QuicKart HydroSurgery System (Includes SpineJet HydroSurgery Console and QuicKart Procedure Cart)	50370		
SpineJet XL Console (Console includes Foot Pedal and Power Cord)	52700		
QuicKart Procedure Cart	50300		
*Blue SpineJet XL Quick Connect must be used with all fusion handpieces. Note, color of Quick Connect must match the color of the handpiece handle.			

All disposable products must be ordered separately.

HydroCision, SpineJet, and The Leader in HydroSurgery are registered trademarks of HydroCision, Inc. ©2004-2006. All rights reserved. Patent No.6,375,635, 6,960,182, 6,923,792, 7,122,017 other US and Foreign Patents Pending. P/N 1000-2378 Rev C 07/09 TG-9904



HydroCision **Customer Care Line** 888.747.4470 Call Toll Free for Product Support

267 Boston Road, Suite 28 North Billerica, MA 01862-2336 Tel. 888.747.4470 Fax 978.600.5058

Email: info@hydrocision.com www.hydrocision.com



The power of water-based technology for better clinical outcomes in minimally invasive spine surgery

HydroSurgery System

SpineJet® XL Advancing fusion disc preparation







Technology advancing fusion disc preparation

Effective

- Accesses entire disc space
- Superior disc prep
- Removes **95%** more posterior contralateral nucleus than conventional devices*

Safe

- 65% fewer instrument passes*
- Unique tip design protects against annular penetration

Efficient

- Aspirates as it resects
- 50% less disc prep time**
- Removes nucleus and cartilage

Reduces risk of:

- Surgical Site Infection (SSI)
- Procedure related complications
- Surgical outcome failure and re-operation
- Exposure to occupational hazards
- Cross contamination

"This single, elegant device accomplishes what it used to take three or four instruments to do, thus empowering us to perform truly MIS fusion procedures that lead to improved clinical outcomes."

> Larry Khoo, M.D. Neurosurgeon The Spine Clinic of Los Angeles Los Angeles, CA

"With minimal passes into and out of the disc space, The SpineJet XL more effectively removes disc material, optimizing the environment for lumbar interbody fusion."

> Mitchell A. Hardenbrook, M.D. Orthopedic Spine Surgeon Boston Spine Group - Boston, MA

*Cadaver studies. Data on file. **User reports. Data on file.

Disc space prepared with conventional tools




Disc space prepared with SpineJet XL



			S
f	•	C	

Measurement*	Standard Tools	SpineJet XL	Benefit
Instrument Passes	102	36	65% Fewer Passes
% Contralateral Tissue Removed	45%	88%	95% More Contralateral Tissue Removed
% Endplate Preparation	70%	86%	23% Greater Surface Area Prepared
% Damaged Endplates	56%	23%	59% Fewer Damaged Endplates
Average Disk Prep Time	17.8 min	9.3 min	50%Time Savings

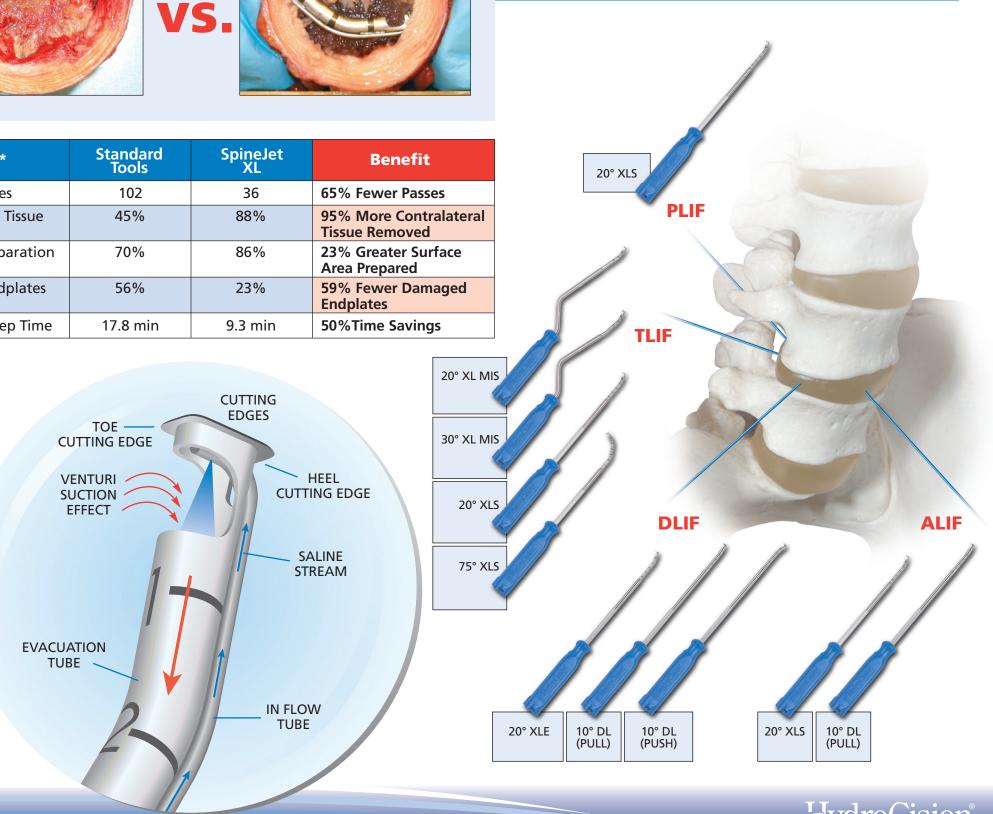

Instrument versatility or all spinal approaches

