

pfmmedical

Quality and Experience

PFO Occlusion System **Nit-Occlud® PFO**

The Nit-Occlud® PFO occluder is a new generation of occlusion systems, leading to optimal performance and effective closure.

Nit-Occlud® PFO is a double-umbrella occluder with dacron membranes in a unique design: The occluder is knitted from a single Nitinol wire. Therefore, it has a low profile and no protruding clamps.

www.pfmmedical.com



Nit-Occlud® PFO is knitted from a single Nitinol wire, making protruding fixation-clamps obsolete. Nit-Occlud® PFO`s unique single-layer distal disc, reduces the metal used in the left-atrium (LA) by 50%, decreasing the thrombo-embolic risk.

Benefits

Patented single-wire knit

Due to the single-wire design, there is no need for protruding clamps which are usually used. This results in a very low profile.

Special design

The occluder consists of a double-layer right atrial disc and a single-layer left atrial disc. The thrombo-embolic risk is lowered by reducing the amount of Nitinol being used in the LA (50%).

Accelerated endothelialisation

A polyester membrane facing the LA promotes accelerated endothelialisation

Various sizes

The Nit-Occlud® PFO is available in diameters of 20 mm, 26 mm and 30 mm.


Easy to use

The easy-to-use application system was developed specifically for Nit-Occlud® PFO procedures. The device is premounted, retrievable and easy to position and to deploy.

Details

- ▶ Single-wire knit
- ▶ No need for protruding clamps to fix loose wire ends
- ▶ Very low-profile
- ▶ Concave single-layer left atrial disc
- ▶ 50% reduction of metal in the left atrium
- ▶ Good adaption to the septal wall
- ▶ Soft and atraumatic rim
- ▶ LA polyester face accelerating endothelialisation
- ▶ Pre-mounted 9 F and 10 F systems
- ▶ Easy to reposition and to retrieve
- ▶ Simple to release
- ▶ Good radiopacity
- ▶ MR conditional

Knowledge

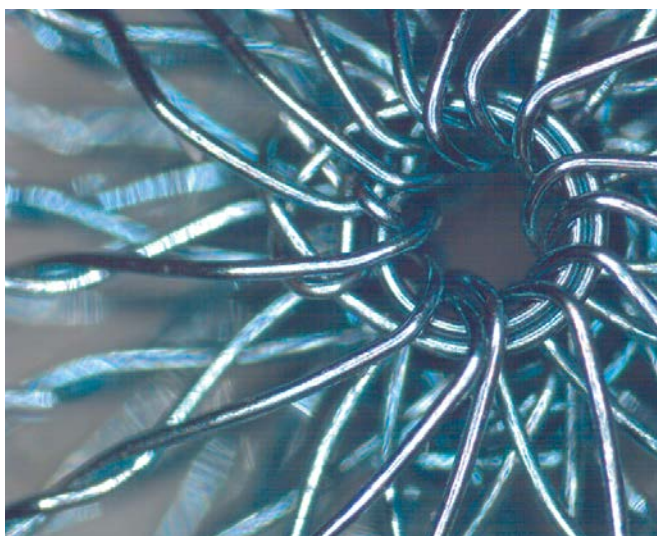

PFO occluders of competitors are made from several wires, with the need to fix the loose ends with clamps.


Nit-Occlud® PFO however is knitted out of a single Nitinol wire, making these fixations obsolete.

Knowledge



The distal part of Nit-Occlud® PFO is designed as a concave single-layer disc. This reduces the amount of metal by half in the left atrium, thus providing best patients benefit by reducing the thrombo-embolic risk.

Function

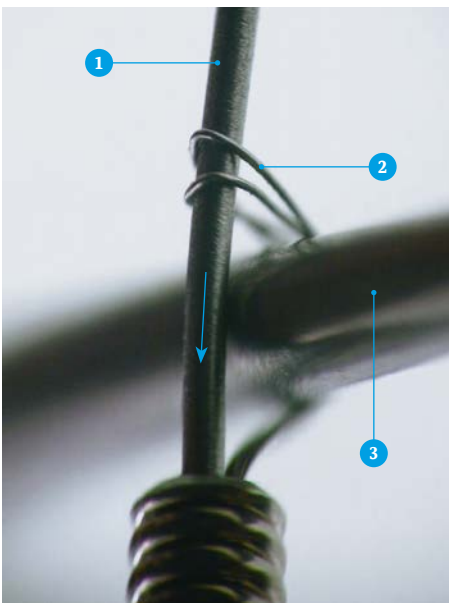

The implant is connected to the pusher via retaining wires. A locking wire fixes this connection. To release the implant, the locking wire has to be pulled through the loop of the retaining wires.

- 1 Locking wire
- 2 Retaining wire
- 3 Part of the implant

Scan the QR code or follow the link and see how the implant works

<http://www.pfmmedical.com/clip-pfo>



Portfolio

Nit-Occlud® Implantation Sheath
PFO Implantation Accessory

Nit-Occlud® Lê VSD
VSD Occlusion Device (coil type)

Multi-Snare®
Grasping device

Nit-Occlud® PDA
PDA Occlusion Device (coil type)

Nit-Occlud® PFO
PFO Occlusion Device (umbrella type)

Ordering Information



Nit-Occlud® PFO › Permanent implant for the closure of Persistent Foramen Ovale (PFO)

› *Content: Implant, loader, delivery system*

REF	Diameter	Recommended Implantation Sheath	Recommended Implantation Sheath Length
182020V1	20 mm	9 F	max. 90 cm
182026V1	26 mm	9 F	max. 90 cm
182030V1	30 mm	10 F	max. 90 cm

Contact

Should you have any questions our Customer Solutions Team will be glad to assist you.

service@pfmmedical.com

+49 (0)2236 9641-221

+49 (0)2236 9641-20

pfm medical ag
Wankelstraße 60
50996 Köln, Germany

Certified according to
DIN EN ISO 13485

Products are subject to alterations
Nit-Occlud®: registered trademark
in Germany

Manufacturer:

pfm medical mepro gmbh, Am Söterberg 4, 66620 Nonnweiler-Otzenhausen, Germany, CE 0124

www.pfmmedical.com - Your information source for the products of pfm medical ag.

PB5005EN/2019-12-04

