


Z-5™

Atrioseptostomy Catheter

**Radiopaque Catheter
Body and Image Band**
for clear visualization and
positioning under fluoroscopy

**35 Degree Tip
Angulation**
facilitates passage into
left atrium

**Non-compliant,
Latex-free balloon**
maintains its shape
during pullback

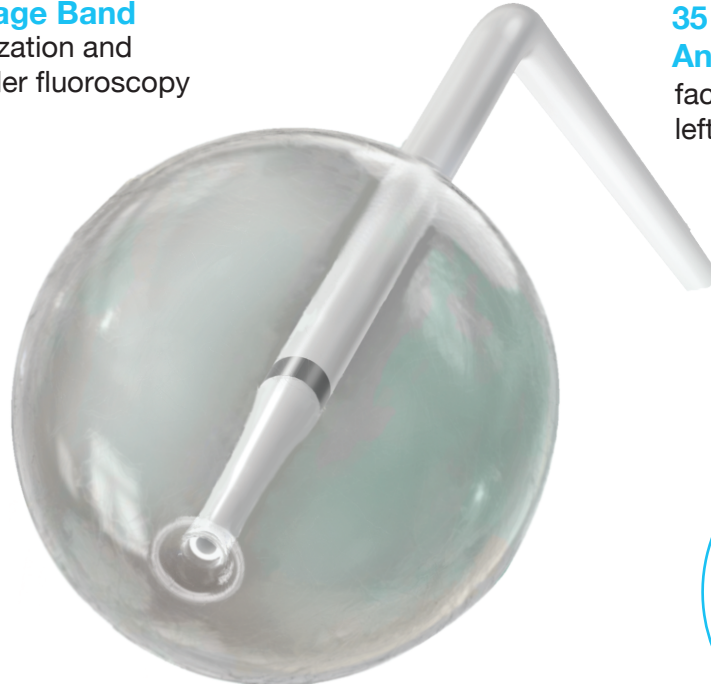

Inner Catheter Lumen
for utilization of a guidewire

Available in
9.5 mm and
13.5 mm
diameters

OVER TWO DECADES OF PROVEN PERFORMANCE

The Z-5™ has been on the market since 1996 and has benefited thousands of patients globally. The Z-5™ is engineered for maximum steering and tracking. The dual-lumen shaft provides pushability coupled with exceptional pull strength. The 9.5mm is primarily for infants less than 2kg.


Indications for use:

Recommended for balloon atrioseptostomy, an accepted technique in most pediatric cardiology centers for the palliation of several congenital cardiac defects.

Refer to the Z-5™ Instructions for Use for full indication statement, relevant warnings, precautions, complications, and contraindications.

Z-5™

Atrioseptostomy Catheter

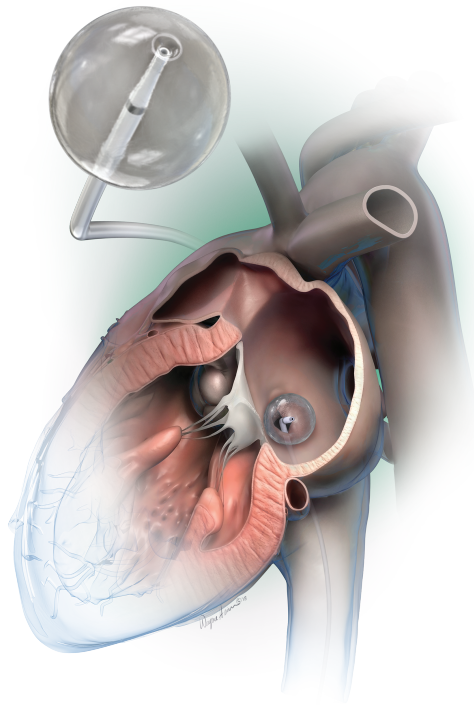

Micro-thin, non-compliant balloon for low deflated profile
Reduced balloon size for neonates with a small left atrium

Z-5™ Specifications

Catalog No.	Balloon Diameter (MM)	Balloon Length (CM)	Introducer Size (FR)	Shaft Size(FR)	Usable Length (CM)	Guidewire (Inches)	Maximum Volume (CC)
SPT002	9.5	.95	5	4	50	.014	1
SPT003	13.5	1.35	6	5	50	.021	2

- Rx only
- Z-5™ is a trademark of NuMED, Inc.

The Z-5™ is available worldwide. Please contact us or your local distributor to place an order.

This device is subject to individual country regulations regarding the importation and sale of this product



G. van Wageningen B.V.
European Representative
info@heartmedical.nl
Hallenweg 40
5683 CT Best
The Netherlands
Tel: +31-499-377388
Fax: +31-499-377456

NuMED | 2880 Main Street | Hopkinton, NY 12965 USA
Tel: 315.328.4491 | Fax: 315.328.4941 | www.numedforchildren.com

CE 1639